



lymphology measuring for custom-made flat-knit lower extremity

effective October 1, 2017

Notes on taking measurements for custom-made flat-knit stockings

• Flat knit compression garments should not be measured until the best possible decongestion effort has been achieved. Circumference () and length () measurements are taken in a distal to proximal sequence (starting at the toes and working up the leg towards the waist).

• It is essential to mark measuring points on the leg so that the circumference and length measurements are taken at the same point.

• The amount of tension applied must be determined by a person experienced in bandaging techniques and depends on edema and tissue characteristics, sensitivity to pain, fibrosis, location of the edema, mobility, and strength of the patient.


• Generally flat-knitted garments are measured under tension with the exception of cY. According to the measuring point we distinguish:

Skin measurements = loose; Measuring to the first resistance = measuring undertension without force; Full Tension = absolutely compression of the tissue.

Knee high

All length measurements are contoured along the medial (inner) side of the leg

Foot measurement points (Y, A, A1, Z)

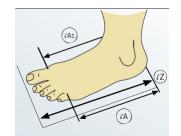

A (for straight and oblique border)

Circumference () Y

Measurements taken over the heel and arch in dorsal flexion with tendons taut. Take physical dimensions without tension. Take soft tissue projections into account.

Length()A

Distance from heel to base of little toe (length of outer side of foot), including the joint.

Length()A1

Distance from heel to base of big toe (medial length), including the joint.


Length()Z

Closed toe: distance from heel to end of longest toe. (Not needed for open toe garments, except open toe seamed toe caps. Z needed for ALL seamed toe cap garments.)

Circumference () A

Apply tape measure around the metacarpophalangeal joints (first joints of toe).

Measuring point B: just above the malleolus (ankle bone)

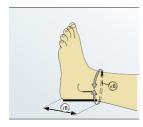


Length () B From sole of foot to above (proximal to)

the ankle .

Circumference () B


Measure circumference, depending on indication, tightly if necessary.

Measuring point C: at the largest calf circumference Measuring point D: at the head of the fibula

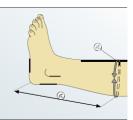

Length () C From sole of foot to largest calf circumference.

Circumference () C Measure circumference, depending on indication, tightly if necessary.

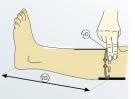

Length()D

From the sole of the foot to the head of the fibula (approx. two finger widths below kneecap).

Circumference () D

Measure circumference, depending on indication, tightly if necessary. Exception: AD stocking = Body dimension.

Length () B1 From sole of f

Measuring point B1: at the base of the calf


From sole of foot to base of the calf.

Circumference () B1 Measure circumference, depending on indication, tightly if necessary.


Thigh high

All length measurements are contoured along the medial (inner) side of the leg except (I)K2 which is contoured on the posterior (back) side of the leg.

Measuring point E: at the middle of the kneecap

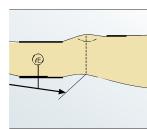


Length () E From sole of foot to middle of kneecap.

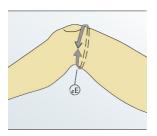

Circumference () E

In general, measurements are not made tightly when the leg is bent (angle 30° - 45°), take circumference measurement one finger width proximal (above) to the back of the knee, and two finger widths proximal (above) to the knee cap. Do not place tape in the back of the knee or folds of skin.

Measuring point F: at the middle of the thigh

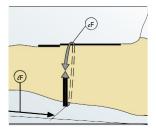

Length () F From sole of foot to middle of thigh.

Circumference () F Depending on the indication, the circumference is generally measured fairly tight.



Measuring point G

Length ()G

AG stocking dimension: sole of the foot to largest point (inner aspect) below the transverse buttock fold (including top band).



Circumference () G


Depending on the indication, the circumference is generally measured fairly tight at the widest part of upper thigh (not always the top of the thigh.)

Measuring point K1: at the pubic bone



Length () K1 From sole of foot to the pubic bone .

Measuring point K2: at the gluteal fold

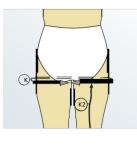

Length () K2 From sole of foot to gluteal fold.

Circumference () K Measure around both legs at the gluteal fold.

Obligatory measurement for medi!

Waist high

All length measurements are contoured along the body except "H" and "T" which are straight lengths.

Measuring point H: widest point of the pelvis




Length () H From sole of foot to widest point of pelvis.



Circumference () H Circumference of pelvis at the widest point, or buttocks, measure tightly depending on indication.

Measuring point T: at the waist

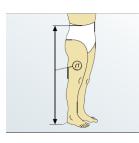

Length () T From the sole of the foot to the top border of the garment (patient-specific).




Circumference () T Circumference of torso, measured under moderate tension.

Measuring points for front () TK1 and back () TK2 top border of the garment



Length () TK1


Front length TK1: measure vertically next to body, at the level of the pubic bone to the top border of the garment.

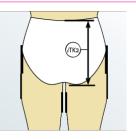

Length () TK2

Back length TK2 measure vertically near the body, at the level of the transverse gluteal fold to the border of the garment.

Measuring a suspensory

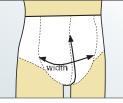

Length & width


Measure the width and length of the suspensory, as marked in Figures A and B - this measurement can be taken over the underwear.



Fig. B

LOWER EXTREMITY MEASURING & ORDER FORM: CUSTOM-MADE FLAT-KNIT



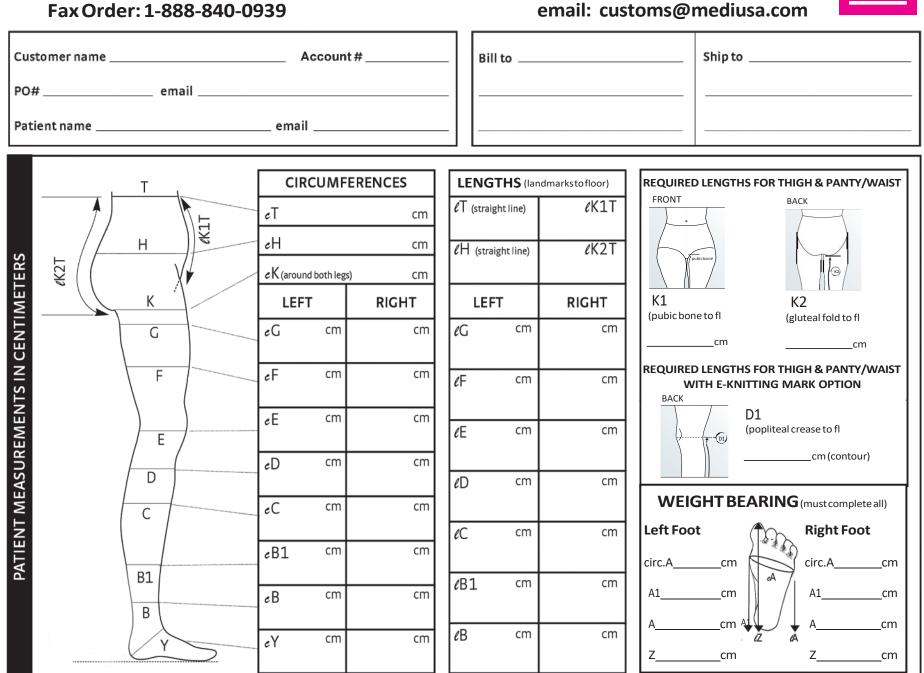

LOWER EXTREMITY MEASURING & ORDER FORM: CUSTOM-MADE FLAT-KNIT

Fax order: 1-888-840-0939 email: customs@mediusa.com

ATERIAL	□ mediven® 550 □ mediven® mondi □ mediven® sensoo (ccLII)		OPTIONS	FOOT Open toe Closed toe: Tricot (half compression) *default		OBLIQUE ENDING Standard Straight	□ Steep	
STYLE/MATERIAL	Knee highWaist high	Vaist high		MONDI SILVER	 Netting (no compression; 550 only) S50 ONLY OPTIONS Gluteal shaper (tricot gusset only) Hallux ease (closed toe or seamless toe cap only) 		ADDED LENGTH IN FRONT (include <code>/K1</code> and <code>/K2</code> measurements below) Needed for all thigh high stockings GUSSET: SUSPENSORY:	
	Men's leotard	Maternity panty * (slight compression panel)		-	PTIONS DY knitting mark at ankle		Tricot (standard)Tricot	
COMPRESSION	PANTY SECTION CCLI: 18-21mmHg CCLIII: 34-46mmHg	CCLII: □ 23-32mmHg CCLIV: □ 49-60mmHg**		E knitting mark at knee SILICONETOPBAND None Profile (honeycomb pattern, 5cm wide only)		NettingNetting CompressiveCompressive Widthcm	Compressive	
	LEFT LEG CCLI: 18-21mmHg CCLIII: 34-46mmHg	CCLII: 23-32mmHg CCLIV: 49-60mmHg **		 Prome (noncycomb pattern, sem wide only) Beaded: Small (2.5cm wide) Wide (5cm wide) Peony (beaded, 5cm wide only) 			Lengthcm LYMPHPADS* Removeable (standard)	
	RIGHT LEG CCLI: 18-21mmHg CCLIII: 34-46mmHg	CCLII: 23-32mmHg CCLIV: 49-60mmHg**		POSITION	SILICONE PIECE	ANTI-SLIP SILICONE DOTS 10 X 4.5 cm	Width (cm) Length (cm) OR Silk Lining Width (cm) Length (cm)	-
QTY	Left	RightPair		Anterior/vertically over "E"	🗖 8 X 5 cm	☐ 10 X 4.5 cm	Location	
COLOR	Standard Colors Caramel Black Cashmere* Sand	Trend Colors* Anthracite Magenta Aqua Mint-Green Blue-Jeans Moss-Green		Posterior/rear over seam Sole of the foot	□ ^{8 X 5 cm} N/A	□ ^{10 X 4.5 cm}	ZIPPER*	or 🗆 Medial 📄 Lateral
DESIGN	mediven 550 Fashion E Mustpick1color&1pattern Colors Patterns □ Berry □ Animal □ Brown □ Crosses	pattern Design Elements ns □ nal □ Ribs		WAISTBAND □ Adjustable (Standard) □ Velcro □ Silicone beaded (5cm wide) □ Knitted endir LEVAMED SILICONE ANKLE PADS* Add 1 cm to cY measurement per pad. Left: □Inner □ Outer □ Permanent □ Removeal			Brilliance	
	Brown Crosses Gray Ornaments					nanent 🛛 Removeable nanent 🖾 Removeable	L Romance	ROMANCE

*Not available in mediven sensoo. Trend colors and fashion elements require an extra 5 days for production. **mediven 550 only

NOTES:

me